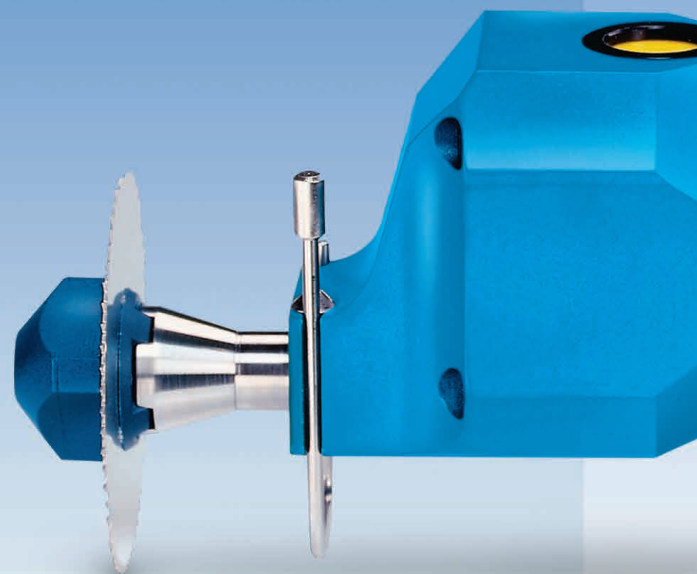




SWISS
MADE

OSCHMED

+ SWISS MEDICAL COMPANY



THE COMPANY OSCIMED SA

Oscimed was founded in 2004. It is based in the swiss Jura, famous area well-known for its watch industry and high technologies sector.

We focus on offering premium medical products for better daily living, with a high quality service, for the satisfaction of our customers.

We pride ourselves on our following qualities:

- Our products carry an official "Made in Switzerland" certification. The Oscimed 2000 and PSV II saws operate with a cutting-edge swiss-made motor and the Ergo II is fitted with a high-quality german motor for optimal performance.
- Availability and flexibility, being a close-knit family company with a focus on trust we aim to work closely with each of our partners in order to find simple, quick and efficient solutions to any problem.

We believe that an integrated system protects the manufacturer while an open system protects the consumer. We therefore provide an open system for consumables and maintenance, which lets you easily swap out blades without any tool and perform maintenance and repairs on-site, additionally if you encounter bigger issues you can send the saw back to Oscimed for quick processing.

Our cast cutter saws are among the lightest and lowest noise options on the market, making them easier to use and more comfortable for your ears.

We are proud to work with valued partners on every continent, several of whom have been with us for nearly three decades, some have had the same Oscimed saw for close to 20 years.

As proof of the unmatched reliability of our products, **a 3-year warranty applies to all our saws !**

Contact us for a quote or any inquiries you might have !



OSCIMED SA – Fritz-Courvoisier 105 – CH-2300 La Chaux-de-Fonds – Tel. +41(0)32 926 63 36 – Fax +41(0)32 926 63 50 – info@oscimedsa.com

TESTIMONIALS

"Oscimed has been a trusted partner of Atlas International for more than a decade, consistently delivering premium Swiss-made cast saws to the U.S. market. Their products are highly regarded for superior quality and reliability, regularly outperforming major competitors and earning strong recognition from hospitals and medical professionals."

Atlas International, United States of America

www.atlasprosales.com



"Your products consistently demonstrate outstanding quality and reliability, with an exceptionally low defect rate. Working with your team has also been a seamless experience thanks to your efficient communication and excellent support."

Muranaka Medical Instruments Co., Ltd,
Japan



www.muranaka.co.jp

"At MEDICIAL, we truly appreciate our wonderful collaboration with OSCIMED and the strong relationship we have built over time based on trust, professionalism, and mutual respect. We are especially pleased with the quality and reliability of OSCIMED's caster saws, which have earned the confidence of both our team and our customers. We sincerely look forward to continuing and growing this successful partnership in the years ahead."

MEDICIAL, Greece



www.medicial.gr

"Medezine Ltd and Oscimed have built a strong and successful partnership over the past 15 years, founded on shared excellence within the casting industry. The outstanding quality and reliability of the products, combined with responsive service and close collaboration, have helped establish a loyal client base across the UK."

Medezine Ltd, England



www.medezine.co.uk

"Our collaboration with Oscimed has always been efficient, professional, and built on excellent communication. Their team is highly supportive, attentive to customer needs, and consistently delivers reliable service and high-quality products. We truly value the strong partnership and look forward to continuing our successful cooperation in the future."

Mida Tecnologia Medica S.P.A., Italy



www.midatecnologiamedica.it

OSCIMED 2000

Cast
cutter
saw



LIGHTWEIGHT AND LOW NOISE, PERFECT FOR CHILDREN PLASTER ROOMS OR HAND SURGEONS

In addition to special saw blades with a noise-suppressing central core, the OSCIMED 2000 cast cutter saw is equipped with a state-of-the-art brushless and fanless SWISS MOTOR for a precise use and one of the lowest noise level among cast cutter saws on the market.

Light, nimble and practical with its power button, this model provides great accuracy and handling.


With a modern shape for single hand use, this ergonomic design combined to SWISS MOTOR technology has proved itself over the years.

Rated power	100 Watts
Length of the saw body	210 mm
Effective diameter	46 mm
Weight without power supply	760 gr.
Speed of motor rotation	13'200 rev/min
Noise level	62 dBA
Length of the power supply cable	1,00 m
Length of the 24 VDC cable	3,50 m

OSCIMED 2000 Saw:

- with European plug: ref. OSC 20101
- with British plug: ref. OSC 20101-UK
- 120V/60 Hz Saw with U.S. plug: ref. OSC 20101-120

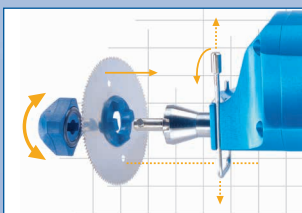3 Years Warranty




Adjustable guide and depth gauge. ref. OSC 203



Usable with one hand thanks to its very short and lightweight shape.



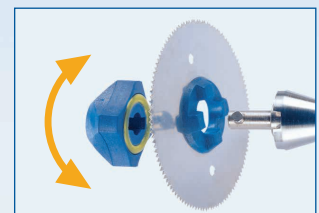Very low noise level thanks to a patented anti-noise technology.



Blade change without tools.

OSCIMED ERGO II

Cast
cutter
saw



OSCIMED
ERGO II



ERGONOMIC, COMPACT, LIGHTWEIGHT AND EASY TO USE

OSCIMED ERGO II is a compact, lightweight plaster saw ergonomically shaped to enable precision work without any effort, even with only one hand. OSCIMED ERGO II is robustly built and equipped with a cooling fan for intensive use without overheating. The device is powered by a safe low voltage of 24 volts and complies with the latest medical standards.

OSCIMED ERGO II is a practical tool for daily use, it allows for quick, accurate and efficient work. It also accommodates a large variety of blades and allows swapping them out with no tools required.

Rated power	180 Watts
Length of the saw body	230 mm
Effective diameter	58 mm
Weight without power supply	860 gr.
Speed of motor rotation	14'500 rev/min
Noise level	67 dBA
Length of the power supply cable	1,10 m
Length of the 24 VDC cable	3,00 m

OSCIMED ERGO II Saw:

- with European plug: ref. OSC 340
- with British plug: ref. OSC 340-UK
- 120V/60 Hz Saw with U.S. plug: ref. OSC 340-US

3 Years Warranty

4



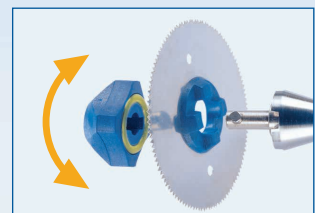
Cooling fan for use without overheating.



High-performance, easy-to-adapt dust vacuuming system.



Usable with one hand thanks to its very short and lightweight shape.



Blade change without tools.

OSCIMED DUST EXTRACTION SYSTEM

4 levels of filtration, adjustable power and automatic start

Vacuum cleaner with adjustable suction power, mounted on four wheels. Its robust and solid construction also provides a platform for the saw.

The HEPA-compliant filters provide a standard filtration of 99.97% of particles of 0.3 microns. A power outlet allows the connection of the saw. The vacuum cleaner starts automatically upon activation of the saw and stops a few seconds after the saw stops.

Power	Adjustable from 300 to 1000 Watts
Noise level	60-70 dB
Suction	1800 mm water column, 40 litres/sec.
Length of the suction hose	2.5 meters.
Capacity of the stainless steel tank	20 litres.
Weight	12 kg.



Vacuum cleaner for OSCIMED ERGO II and OSCIMED 2000

OSCIMED vacuum cleaner with European plug: ref. OSC 206-V90

OSCIMED vacuum cleaner with British plug: ref. OSC 206-V90-UK

OSCIMED 120 V/60 Hz vacuum cleaner with U.S. plug: ref. OSC 206-V90-US

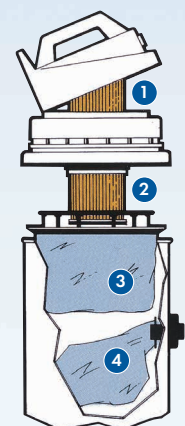


Storage compartment. Handles make it easy to move.

Filtration system with 4 levels

Standard HEPA: 99.97% of particles of 0.3 microns

- 1 HEPA filter ref. OSC 208-3
- 2 HEPA filter ref. OSC 208-2
- 3 Polyester filter ref. OSC 208-1
- 4 Paper dust bag ref. OSC 207





3 Years Warranty

SWISS MOTOR INSIDE, BEST TECHNOLOGIES FOR EASY USE

Thanks to its small and powerful SWISS MOTOR, OSCIMED PSV II is a very well balanced and miniaturized cast cutter saw. Its light weight allows single-handed use even for nurses. Its ergonomic design includes a grip area with reduced circumference. The motor is mounted on rubber dampers to reduce noise and vibrations. The cooling ventilator permits intensive use without overheating. Accessories such as blades and vacuum nozzles can be changed rapidly, without tools. The device operates on a safe low voltage and meets the latest medical standards. This saw is delivered with our special anti-noise blade BS-65-LN.



Anti-noise blade BS-65-LN



Installation guide

AVAILABLE IN TWO VERSIONS, WITH OR WITHOUT DUST EXTRACTION

The PSV II without air duct for use without vacuum cleaner is a VERY smart saw with one of the smallest handle circumference. This model is the new reference on the market. The PSV II is perfect for children plaster rooms or hand surgeons. The PSV II with air duct is made for use with a vacuum cleaner. Our engineers have designed this model to minimize handle size and greatly improve freedom of movement.

	with air duct	without air duct
Rated power	180 Watts	180 Watts
Length of the saw body	210 mm	210 mm
Handle size	49x48 mm	48x34 mm
Weight without power supply	755 gr.	640 gr.
Speed of motor rotation	14'100 rev/min	14'100 rev/min
Noise level	63 dBA	63 dBA
Length of the power supply cable	1,80 m	1,80 m
Length of the 24 VDC cable	3,00 m	3,00 m

OSCIMED PSV II Saw with air duct

with European plug: ref. OSC 400V

with British plug: ref. OSC 400V-UK

120V/60 Hz Saw with U.S. plug: ref. OSC 400V-US

OSCIMED PSV II Saw without air duct

with European plug: ref. OSC 400

with British plug: ref. OSC 400-UK

120V/60 Hz Saw with U.S. plug: ref. OSC 400-US



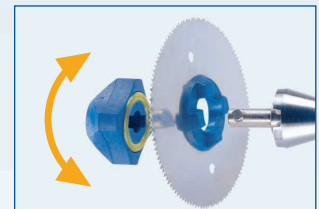
PSV II with air duct.



PSV II without air duct.



Usable with one hand thanks to its very short and lightweight shape.



Blade change without tools.

OSCIMED DUST EXTRACTION SYSTEM

4 levels of filtration, adjustable power and automatic start

Vacuum cleaner with adjustable suction power, mounted on four wheels. Its robust and solid construction also provides a platform for the saw. A power outlet allows the connection of the saw. The vacuum cleaner starts automatically upon activation of the saw and stops a few seconds after the saw stops.

Power	Adjustable from 300 to 1000 Watts
Noise level	60-70 dB
Suction	1800 mm water column, 40 litres/sec.
Length of the suction hose	2.5 meters.
Capacity of the stainless steel tank	20 litres.
Weight	13 kg.



The saw and the vacuum hose can rotate independently. No electric cable attached to the hose or laying on the floor.



Storage compartment. Handles make it easy to move.



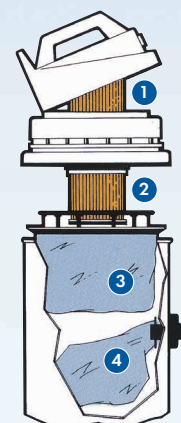
Vacuum cleaner for OSCIMED PSV II

OSCIMED vacuum cleaner with European plug: ref. OSC 206-V00B
 OSCIMED vacuum cleaner with British plug: ref. OSC 206-V00B-UK
 OSCIMED 120V/60 Hz vacuum cleaner with U.S. plug: ref. OSC 206-V00B-US

Filtration system with 4 levels

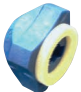







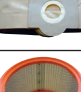



Standard HEPA:
 99.97% of particles of 0.3 microns

- 1 HEPA filter ref. OSC 208-3
- 2 HEPA filter ref. OSC 208-2
- 3 Polyester filter ref. OSC 208-1
- 4 Dust bag ref. OSC 207



SAW BLADES AND ACCESSORIES

Saw Blades	Ø 77 mm	Ø 65 mm	Ø 65 mm	Ø 65 mm	Ø 50 mm	Ø 50 mm
Stainless steel saw blade with white core for plaster of Paris casts.	–	WP-65	–	WP-65-2	WP-50	–
Hardened steel saw blade with blue core for synthetic casts.	BS-77-2	BS-65	BS-65-LN	–	BS-50	BS-50-2
Hardened steel saw blade with blue core with Titan (TiN) processing for synthetic casts.	–	BS-65-GT	BS-65-LN-GT	–	BS-50-GT	–
Hardened steel saw blade with blue core with non-stick PTFE processing for synthetic casts.	–	BS-65-TF	–	–	–	–

Accessories		Cast cutter saw					Vac. cleaner
Description	Reference	Osc 2000 OSC 20101	Ergo II OSC 340	PSV OSC 300	PSV II OSC 400	PSV II with Vacuum OSC 400V	Dust extraction OSC 206-XX
 Blade lock nut	OSC 202	X	X	X	X	X	
 Vacuum nozzle	OSC 240	X	X				
 Vacuum nozzle	OSC 320			X			
 Vacuum nozzle	OSC 420					X	
 Vac. clean. Hose	OSC 272	X	X				
 Vac. clean. Hose	OSC 311			X			
 Vac. clean. Hose	OSC 411					X	
 Polyester filter	OSC 208-1						X
 Dust bag (10x)	OSC 207						X
 HEPA filter	OSC 208-2						X
 HEPA filter	OSC 208-3						X
 Adjustable guide	OSC 203	X					